

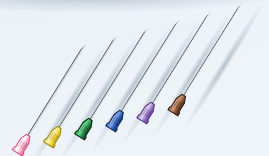
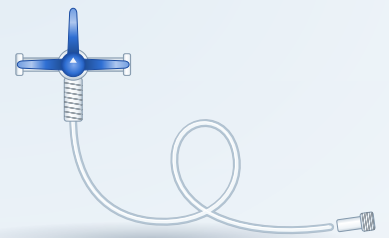
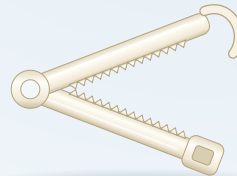
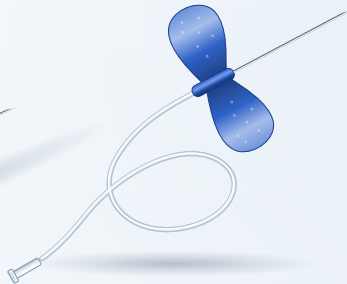
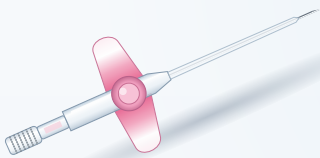
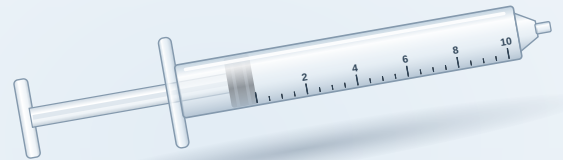


Disposable Medical Devices

Sterile, single-use syringes, needles, infusion sets, catheters and neonatal care devices.

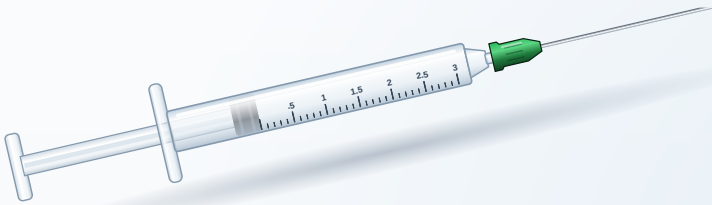

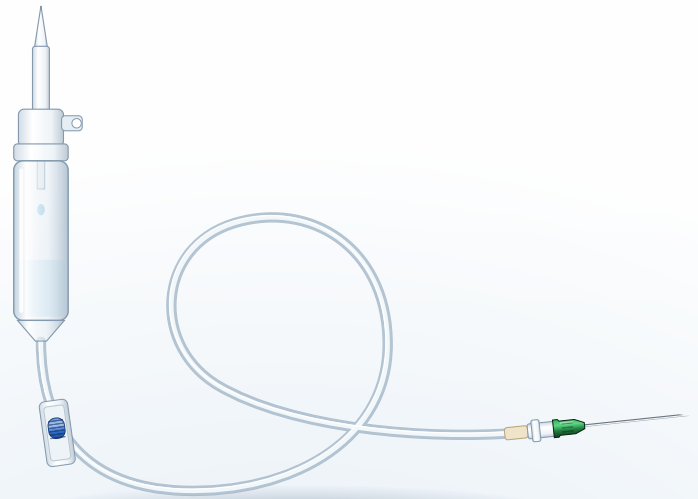

Product range

17 devices in five categories. Sizes shown are the full licensed range; see each product page for details.

Product	Size range	Class	Shelf life	Page
Syringes				
Disposable Hypodermic Syringe with Needle	0.5 ml to 50 ml	B	5 years	3
Disposable Hypodermic Syringe without Needle	0.5 ml to 50 ml	B	5 years	3
Auto-Disable Syringe with Needle	0.5 ml to 50 ml	B	5 years	3
Syringe packing and carton guide	2 ml to 10 ml			4
Needles				
Disposable Hypodermic Needle	16G to 30G	B	5 years	5
Spinal Needle	14G to 26G	B	5 years	5
Blood Collection Needle	20G to 22G	B	3 years	5
Infusion Therapy				
Infusion Set (IV Set), Economy and Premium	21G × 1.5" needle	B	5 years	6
Blood Transfusion Set (BT Set)	18G needle	B	5 years	6
Scalp Vein Set	18G to 25G	B	5 years	6
IV Cannula	14G to 26G	B	5 years	7
3-Way Stop Cock	Line 10 cm to 250 cm	B	5 years	7
Extension Tube	10 cm to 250 cm	B	5 years	7
Tubes & Catheters				
Ryles Tube (Nasogastric Tube)	8 FG to 20 FG	B	5 years	8
Foley Catheter	300 mm to 400 mm	B	5 years	8
Yankauer Suction Set	1/4" and 3/16" bore	B	5 years	8
Neonatal Care				
Infant Feeding Tube	100 mm to 500 mm	B	5 years	9
Umbilical Cord Clamp	All sizes	A	5 years	9

Every product is available ethylene-oxide (EO) sterile. All products except the auto-disable syringe are also available non-sterile. Device class is as per the Medical Devices Rules, 2017(India).



Disposable Hypodermic Syringe with Needle

For transferring solutions and administering drugs. Supplied with the needle fitted, ready for use.

- Crystal-clear barrel with bold, accurate graduations
- Smooth, low-friction plunger movement for precise dosing
- Sharp triple-bevel needle for reduced penetration force
- Manufactured to ISO 7886-1

Sizes	0.5 ml1 ml2 ml3 ml5 ml10 ml20 ml50 ml
Needle	All gauges and lengths
Material	Barrel and plunger: medical-grade polypropylene. Gasket: Santoprene (TPE). Needle: SUS-304 stainless steel
Packing	Ribbon pack or blister pack



Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B

Disposable Hypodermic Syringe without Needle

The same syringe without a needle, for use with the needle of your choice or for transferring solution between containers.

- Universal tip fits all standard needles
- Leak-free gasket seal across the full volume range
- Clear graduations for confident measurement

Sizes	0.5 ml1 ml2 ml3 ml5 ml10 ml20 ml50 ml
Material	Barrel and plunger: medical-grade polypropylene. Gasket: Santoprene (TPE)
Packing	Ribbon pack or blister pack



Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B

Auto-Disable Syringe with Needle

The plunger locks after one full injection, preventing reuse. Made for immunisation and public health programmes.

- Automatic lock engages after a single full stroke
- Prevents reuse and helps curb blood-borne infections
- Supplied sterile, with needle fitted

Sizes	0.5 ml1 ml2 ml3 ml5 ml10 ml20 ml50 ml
Needle	All gauges and lengths
Material	Barrel and plunger: medical-grade polypropylene. Gasket: Santoprene (TPE). Needle: SUS-304 stainless steel
Packing	Individually packed. Bulk and hospital packs available



Sterility	EO sterile only
Shelf life	5 years
Device class	B



Syringe packing

Every syringe is sealed in its own printed pack, then packed in boxes colour-coded by volume so stores teams can pick by sight.



Unit pack

Available as a ribbon pack or blister pack. Each unit pack is printed with:

- Syringe volume
- ISO 7886-1 reference
- Sterilisation method
- Batch number
- Manufacturing and expiry dates
- Licence number

Box colour code and quantities

Single Use Syringe with Needle
STERIX GLOBAL 3
100 pcs | Ribbon Pack
STERILE EO 3 ml/cc

3 ml	Purple box
Per box	100 pcs
Per outer carton	3,000 pcs

Single Use Syringe with Needle
STERIX GLOBAL 5
100 pcs | Ribbon Pack
STERILE EO 5 ml/cc

5 ml	Green box
Per box	100 pcs
Per outer carton	2,000 pcs

Single Use Syringe with Needle
STERIX GLOBAL 10
50 pcs | Ribbon Pack
STERILE EO 10 ml/cc

10 ml	Red box
Per box	50 pcs
Per outer carton	1,000 pcs

2 ml syringes: 3,000 pcs per outer carton. Quantities are for ribbon-packed syringes with needle.

Carton marking

Front panel

Brand, device name, volume, pack type, sterilisation method and QR code

Side panel

Needle size, batch number, manufacturing and expiry dates, licence number and quantity

Private label

Unit packs and cartons printed to your artwork, subject to minimum order quantity



Disposable Hypodermic Needle


For injecting fluids into, or withdrawing fluids from, the body below the surface of the skin.

- Polished, siliconised cannula for smoother insertion
- Colour-coded hub for instant gauge identification
- Available in all lengths across the gauge range

Gauges 16G to 30G (see colour code below)

Material Cannula: SUS-304 stainless steel. Hub and cover: medical-grade polypropylene

Packing Individually packed. Bulk and hospital packs available

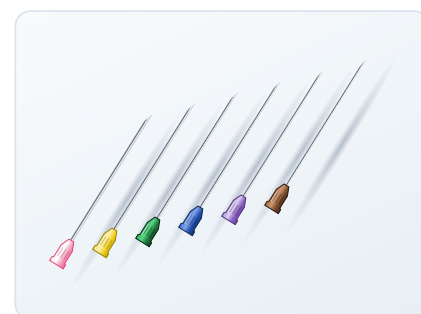

Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B

Hub colour code (ISO 6009)

16G	17G	18G	19G	20G	21G	22G	23G	24G	25G	26G	28G	29G	30G
White	Red-violet	Pink	Cream	Yellow	Deep green	Black	Deep blue	Purple	Orange	Brown	Blue-green	Red	Yellow

Spinal Needle


For injecting local anaesthetic into the subarachnoid space for spinal anaesthesia, or collecting CSF in diagnostic lumbar puncture.

- Precision-ground tip for controlled dural puncture
- Transparent hub shows CSF flashback instantly
- Snug-fitting stylet for smooth insertion

Gauges 14G 16G 17G 18G 20G 22G 24G 26G

Options With or without introducer

Material Cannula: SUS-304 stainless steel. Hub and cover: medical-grade polypropylene

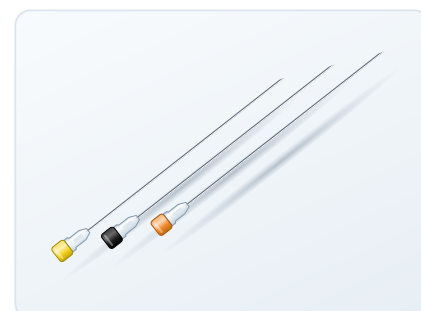

Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B


Blood Collection Needle

For collecting measured volumes of blood into evacuated tubes of various types and capacities.

- Double-ended design for evacuated tube systems
- Self-sealing rubber sleeve prevents leakage between tubes
- Designed for minimal haemolysis and patient discomfort

Gauges 20G 21G 22G

Material Polypropylene, stainless steel, silicone, rubber sleeve

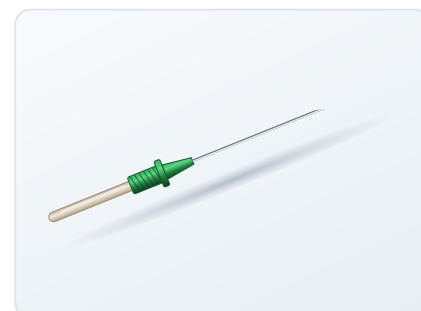

Sterility	EO sterile / Non-sterile
Shelf life	3 years
Device class	B



Infusion Set (IV Set, Vented)

For gravity-feed IV therapy, or with a pressure infuser where an extended fluid path is needed. Use with an IV cannula; the attached needle can also go directly into a vein.

- Transparent drip chamber for accurate flow monitoring
- Precision roller clamp for fine flow control
- Vented spike fits all standard fluid containers


Economy IV Set

Printed pouch

IV Set Premium

Individual blister pack

Needle	21G × 1.5" standard. Other gauges from 16G to 30G on request
Drop rate	20 drops ≈ 1 ml
Material	Polypropylene, ABS, soft PVC, HIPS, latex rubber, SUS-304 stainless steel

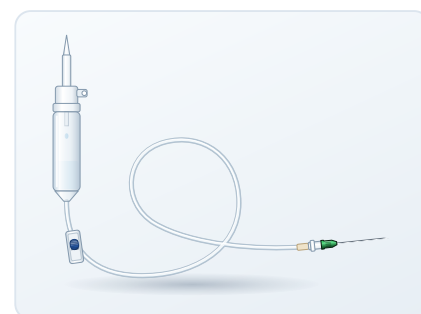


Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B

Blood Transfusion Set (BT Set, Vented)

For delivering measured, regulated blood and blood components to the patient.

- 200-micron cellulose filter chamber traps clots and debris
- Sharp vented spike pierces blood bags cleanly
- Roller regulator for measured, controlled delivery

Needle	18G standard. Other gauges from 14G to 26G on request
Material	Polypropylene, ABS, soft PVC, HIPS, latex rubber, stainless steel, cellulose filter
Packing	Individually packed. Bulk and hospital packs available

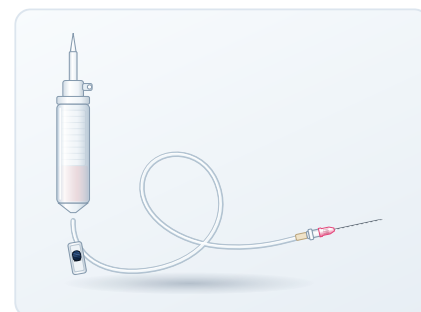


Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B

Scalp Vein Set

Winged butterfly set for routine blood collection and short-term intravenous infusion.

- Soft butterfly wings give a low-angle, stable grip
- Flexible tubing reduces movement at the puncture site
- Suited to paediatric, geriatric and difficult venous access

Gauges	18G 19G 20G 21G 22G 23G 24G 25G
Material	Polypropylene, stainless steel, PVC, ABS, HDPE, POM
Packing	Individually packed. Bulk and hospital packs available

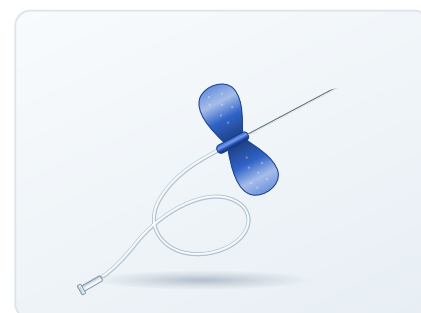

Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B




IV Cannula

Winged cannula with injection port, for infusing fluids, drugs and blood components, and for placing vascular access devices.

- Flexible wings for secure, comfortable fixation
- Integrated injection port for intermittent medication
- Colour-coded by gauge to international convention

Gauges 14G to 26G (see colour code below)

Material Polypropylene, ABS, soft PVC, HIPS, stainless steel

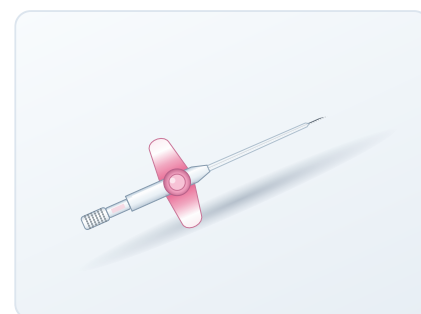

Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B

Colour code (ISO 10555-5)



14G

Orange



16G

Grey



17G

White



18G

Green



20G

Pink



22G

Blue



24G

Yellow



26G

Violet

3-Way Stop Cock


Directional control of fluid flow, with an access port for giving solutions, withdrawing fluid and pressure monitoring.

- Smooth 360° rotating handle with clear flow indication
- Leak-free luer connections on all three ports
- Lipid-resistant version for nutrition therapy

Versions Lipid and non-lipid

Options With or without extension line, 10 cm to 250 cm

Material LDPE, HDPE, polypropylene, PVC, K-resin

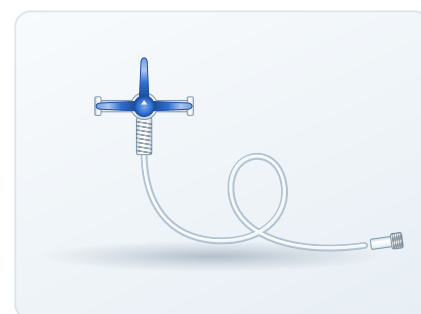

Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B


Extension Tube

Single-use extension line, used as part of a system for infusing fluids and medication.

- Kink-resistant tubing keeps flow uninterrupted
- Secure luer-lock connections at both ends

Length (cm)	10	15	20	25	30	50	80	100	125	150
	175	200	250							

Material LDPE, HDPE, polypropylene, PVC, K-resin

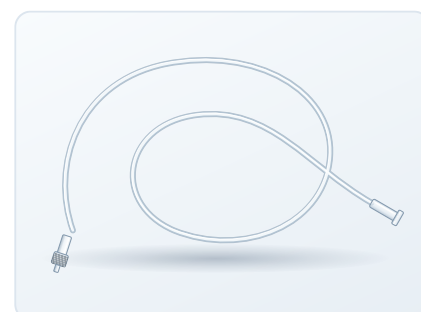

Sterility	EO sterile / Non-sterile
Shelf life	5 years
Device class	B



Ryles Tube (Nasogastric Tube)

Carries food and medicine to the stomach through the nose, for complete or supplementary feeding.


- Radiopaque line for X-ray confirmation of position
- Smooth, rounded closed tip with lateral eyes
- Graduated depth marks for accurate placement

Sizes

[8 FG](#)[10 FG](#)[12 FG](#)[14 FG](#)[16 FG](#)[18 FG](#)[20 FG](#)

Material

PVC and soft PVC with barium sulphate line

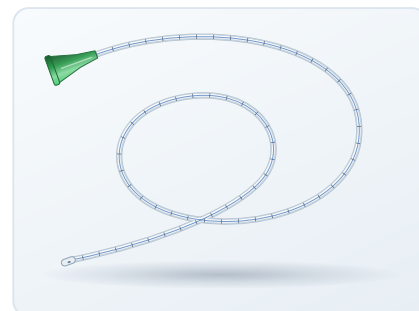

Sterility

EO sterile / Non-sterile

Shelf life

5 years

Device class

B

Foley Catheter

Indwelling urinary catheter, held in the bladder by a balloon inflated with sterile water, for continuous drainage.


- Symmetrical, high-integrity retention balloon
- Smooth surface for comfortable placement
- Colour-coded funnel for size identification

Length

300 mm to 400 mm

Material

Latex rubber, polypropylene

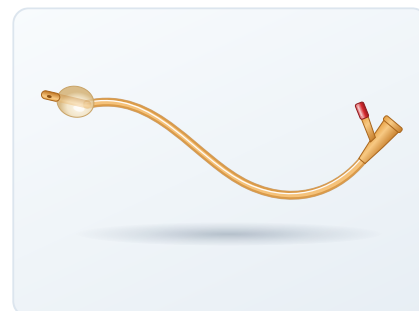

Sterility

EO sterile / Non-sterile

Shelf life

5 years

Device class

B

Yankauer Suction Set

Suction set for ENT procedures, removing excess fluid and debris from the surgical site.

- Ergonomic rigid handle for controlled suction
- Large-bore tubing for rapid fluid removal
- Multiple tip eyes help prevent tissue grab

Bore


[1/4"](#)[3/16"](#)

Tube length

[6 ft](#)[9 ft](#)[10 ft](#)[12 ft](#)

Material

Medical-grade PVC, silicone or polyurethane

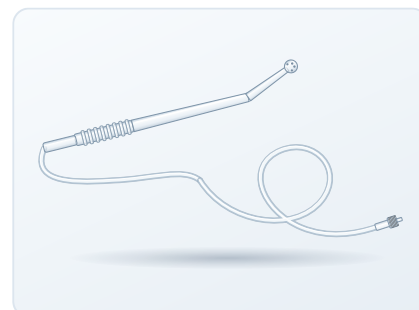

Sterility

EO sterile / Non-sterile

Shelf life

5 years

Device class

B




Infant Feeding Tube

Soft tube placed through the nose or mouth into the stomach, giving feeds and medicines until the baby can take food by mouth.

- Ultra-soft, atraumatic closed tip
- Clearly graduated for precise insertion depth
- Designed for neonatal and NICU feeding

Length	100 mm to 500 mm
Material	Natural rubber, polypropylene, PVC
Packing	Individually packed. Bulk and hospital packs available

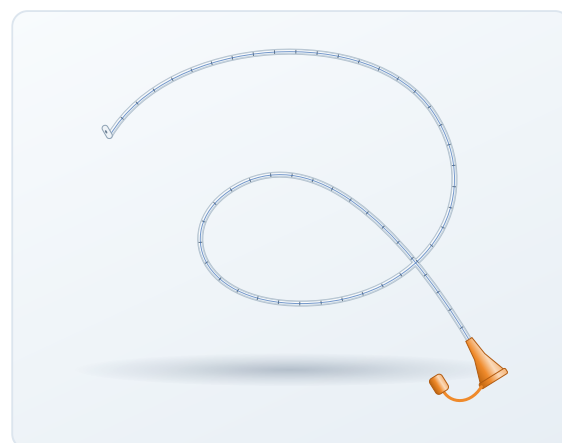


Sterility	E0 sterile / Non-sterile
Shelf life	5 years
Device class	B

Umbilical Cord Clamp

Holds the umbilical cord securely when the cut is made, and stays in place for several days while the remaining cord dries.

- Double-grip jaw with secure snap lock
- Interlocking teeth grip without cutting tissue
- Smooth, rounded edges, gentle on newborn skin

Sizes	All sizes
Material	ABS, soft PVC, polypropylene, HIPS
Packing	Individually packed. Bulk and hospital packs available

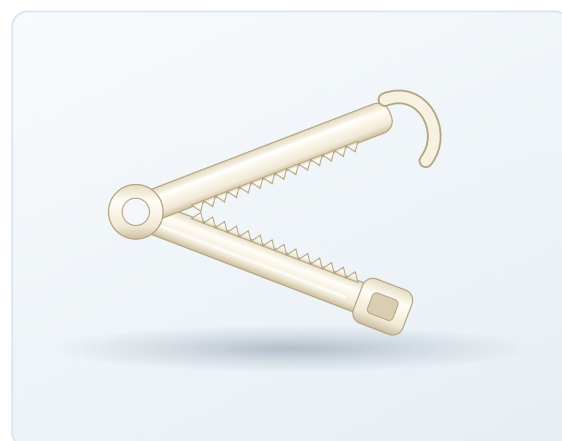

Sterility	E0 sterile / Non-sterile
Shelf life	5 years
Device class	A



Enquiries and samples

Samples and technical documentation are available for every product in this catalogue. To get a quote, send us:

- ✓ Product name
- ✓ Size or gauge
- ✓ Sterile or non-sterile
- ✓ Pack type
- ✓ Quantity
- ✓ Private label, if required

Telephone

+91 84477 05170

Customer care

+91 63884 35686

Manufacturing unit

**Khasra No. 22, Village Palta, Deva Road,
Barabanki, Uttar Pradesh 225301, India**

Registered office

**Flat No. 502, Tower 1B, Sunbreeze 1, BBD
Green, Lucknow, Uttar Pradesh 226028, India**